


Single Use 2-Lumen Sphincterotome V

## CleverCut 2V™

Single Use 3-Lumen Sphincterotome V

## CleverCut 3V™

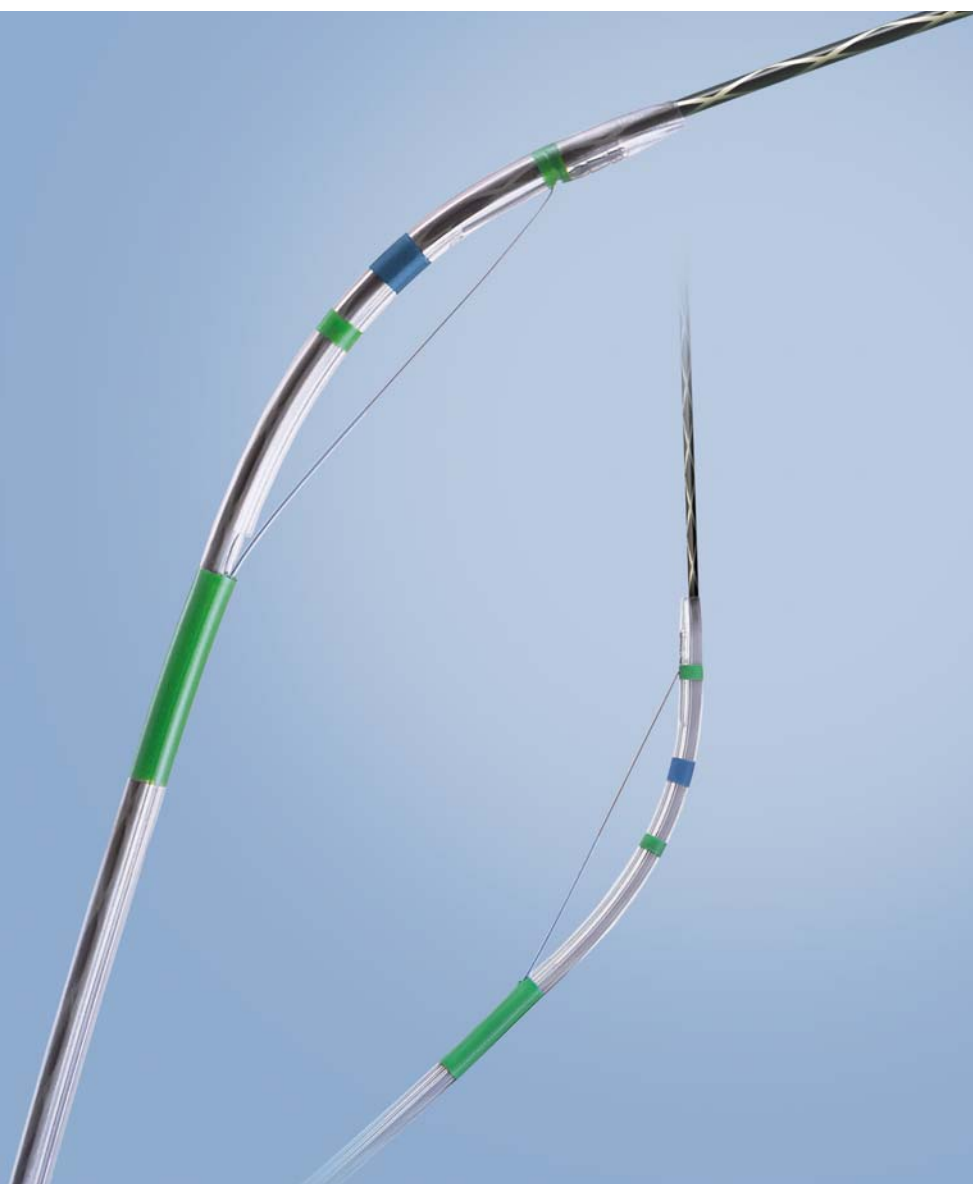


# Enhanced Design for Optimal Cannulation and Minimal Tissue Damage

CleverCut 2 V<sup>™</sup>    CleverCut 3 V<sup>™</sup>

**CleverCut<sup>™</sup> coating enhances safety**

The CleverCut coating on the proximal end of the cutting wire minimizes damage to the surrounding tissue as well as the risk of electrical contact between the wire and endoscope.

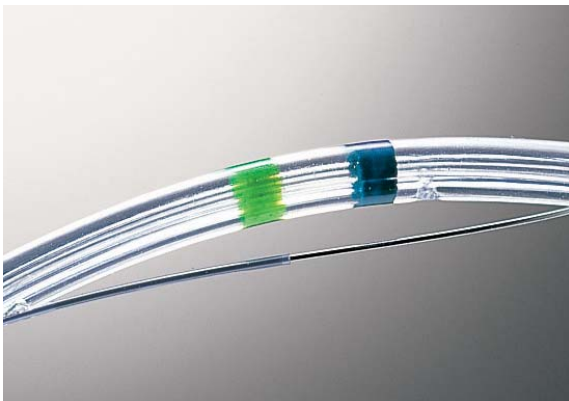

**Radiopaque markings for clear visibility under fluoroscopy**



The radiopaque markers on the CleverCut sphincterotomes provide excellent visibility under fluoroscopy.

**Pre-curved tip to assist with positioning**

The distal end of the CleverCut sphincterotomes is pre-curved to assist with cannulation of the papilla.

**Endoscopic markers assist with positioning**

Color-coded markings are easily visible under endoscopic view and show the best cutting position of the wire.

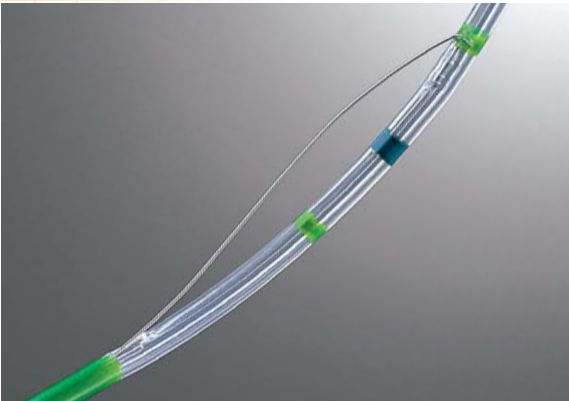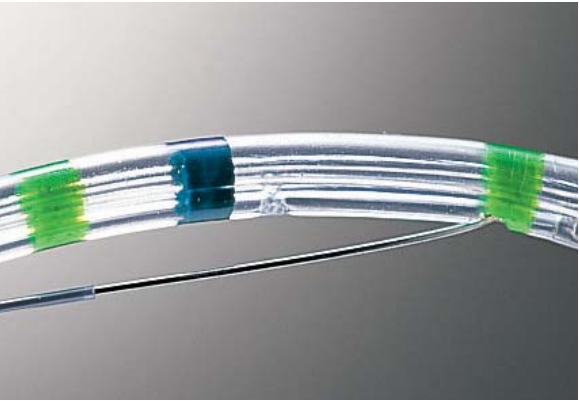

# Single Use 2-Lumen Sphincterotome V

# Single Use 3-Lumen Sphincterotome V



C-Hook



V-Marking




V-Sheath


## Specifications

Model	Lumen type	Distal tip diameter	Maximum insertion portion diameter	Working length	Compatible guidewire	Tip length	Knife length
KD-V411M-0320	3-Lumen type	4.5 Fr	2.5 mm	1700 mm	0.89 mm (0.035 inch)	3 mm	20 mm
KD-V411M-0330						3 mm	30 mm
KD-V411M-0720						7 mm	20 mm
KD-V411M-0725						7 mm	25 mm
KD-V411M-0730						7 mm	30 mm
KD-V411M-1520						15 mm	20 mm
KD-V411M-1530						15 mm	30 mm
KD-V411M-3030						30 mm	30 mm
KD-V431M-0720*		4.0 Fr			0.64 mm (0.025 inch)	7 mm	20 mm
KD-V431M-0730*						7 mm	30 mm
KD-V211M-0720	2-Lumen type	4.5 Fr	2.1 mm	1700 mm	0.89 mm (0.035 inch)	7 mm	20 mm
KD-V211M-0725						7 mm	25 mm
KD-V211M-0730						7 mm	30 mm
KD-V211M-1520						15 mm	20 mm
KD-V211M-1530						15 mm	30 mm
KD-V211M-3030						30 mm	30 mm

\* Tapered tip design for minor or constricted papilla

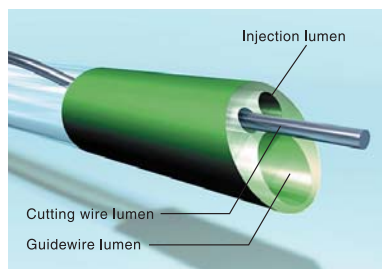

3-Lumen type

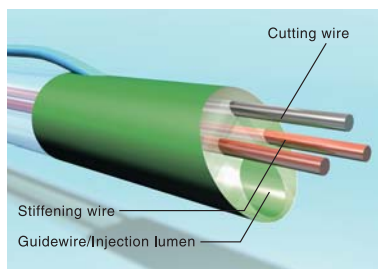

2-Lumen type

Specifications, design and accessories are subject to change without any notice or obligation on the part of the manufacturer.

**OLYMPUS**

OLYMPUS MEDICAL SYSTEMS CORP.

Shinjuku Monolith, 2-3-1 Nishi-Shinjuku, Shinjuku-ku, Tokyo 163-0914, Japan

For a complete listing of  
sales and distribution locations visit:  
[www.olympus.com](http://www.olympus.com)

Printed in Japan F0265E4-022012